




nesplora

 GIUNTI PSYCHOMETRICS

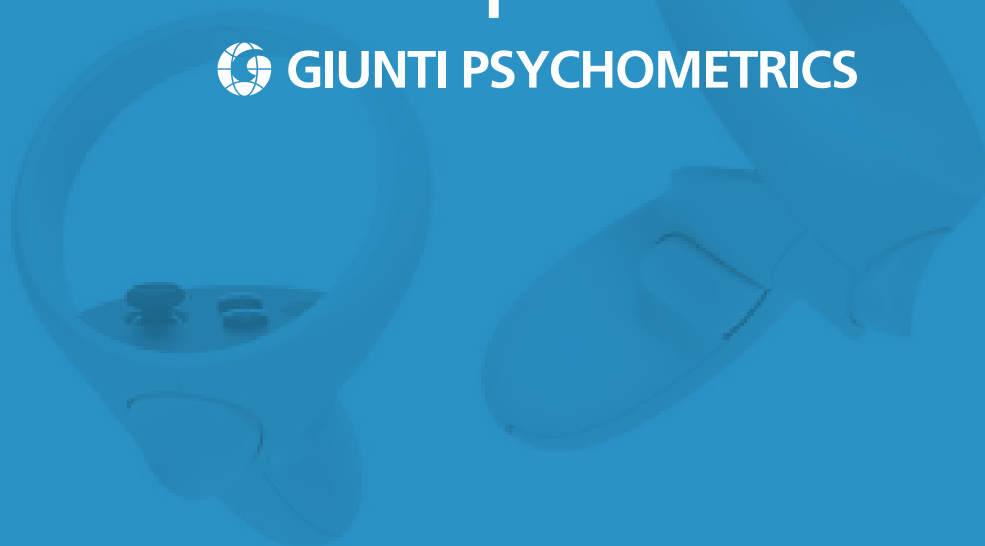

# MINIMUM HARDWARE REQUIREMENTS

In order to use the Nesplora system tests you need to have 4 items:

1. Virtual Reality Device like [Meta Quest 2](#), [Meta Quest 3](#) o [Meta Quest 3S](#)
2. Computer, desktop or laptop
3. Wi-Fi network with an internet connection
4. Wired headset, with cable

NESPLORA DOES NOT SELL THE HARDWARE

## SPECIFICATIONS

### 1 Virtual Reality Device

[Meta Quest 2](#); [Meta Quest 3](#) o [Meta Quest 3S](#)

Available with different capacities (64, 128, 256 ó 512GB) depending on the model. The difference in capacity allows you to install more or less applications on the device, but **any of these are perfectly suitable for any Nesplora test**, which do not use up more than 1 GB on a single device.

#### where can I buy them?

- On the [Official META Manufacturer Website](#)
- At their [Authorised Retailers](#), you can check by country or region



## 2 Computer, either desktop or laptop

With the following characteristics

- **Operating System:** Windows 10 or above, or Mac OS X or above
- **RAM 4GB Memory**
- **Graphics card:** Intel HD Graphics or above

No mobile phone is required to carry out Nesplora System tests.

## 3 Wi-Fi network

With an internet connection

**Virtual reality glasses and computers are linked through a Wi-Fi network.**

When doing the tests and processing the report, both the virtual reality glasses and the computer must be connected to the same Wi-Fi network with an internet connection.

## 4 Wired headset

With cable

**It must cover the ear auricle to sufficiently insulate from any external noise.**

Easily adjustable to adapt to different users.

With leather or similar pads so that they are easier to clean and maintain.

They will be connected to the Virtual Reality device by cable. The type of input, depending on the model is:

- For Meta Quest 2 and Meta Quest 3 → stereo Jack 3,5mm (most common connector for headsets, commonly used with mobile phones and computers)
- For Meta Quest 3S → Type C (the same headset with 3.5 mm stereo jack plug can be used with an adapter that converts from stereo jack to type C.)


**The use of a headset connected by cable to the virtual reality device is mandatory to ensure auditory immersion and stimulus presentation times.**

**Wireless headsets cannot be used for our tests.**



CONTACT  
OUR SPECIALISTS

alobo@nesplora.com  
+34 647 79 75 54  
Online booking [HERE](#)

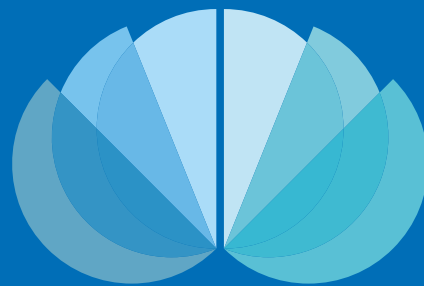

nesplora

PART OF

 **GIUNTI**  
PSYCHOMETRICS